

Safe2 Rotator™ Luer Connector


The connection you can rely on

Easy attach – Firm lock – 360° Rotational flexibility - Easy detach

Good and safe connection of IV tubing is critical to safe and effective IV therapy. IV tubing disconnections can result in serious and irreversible patient injury or even death.

When it comes to stopcocks, apart from being securely connected, they must be convenient to use and perform their primary functions – direct fluids and provide access to the IV set.

With this in mind, Elcam Medical has designed the new Safe2 Rotator™.

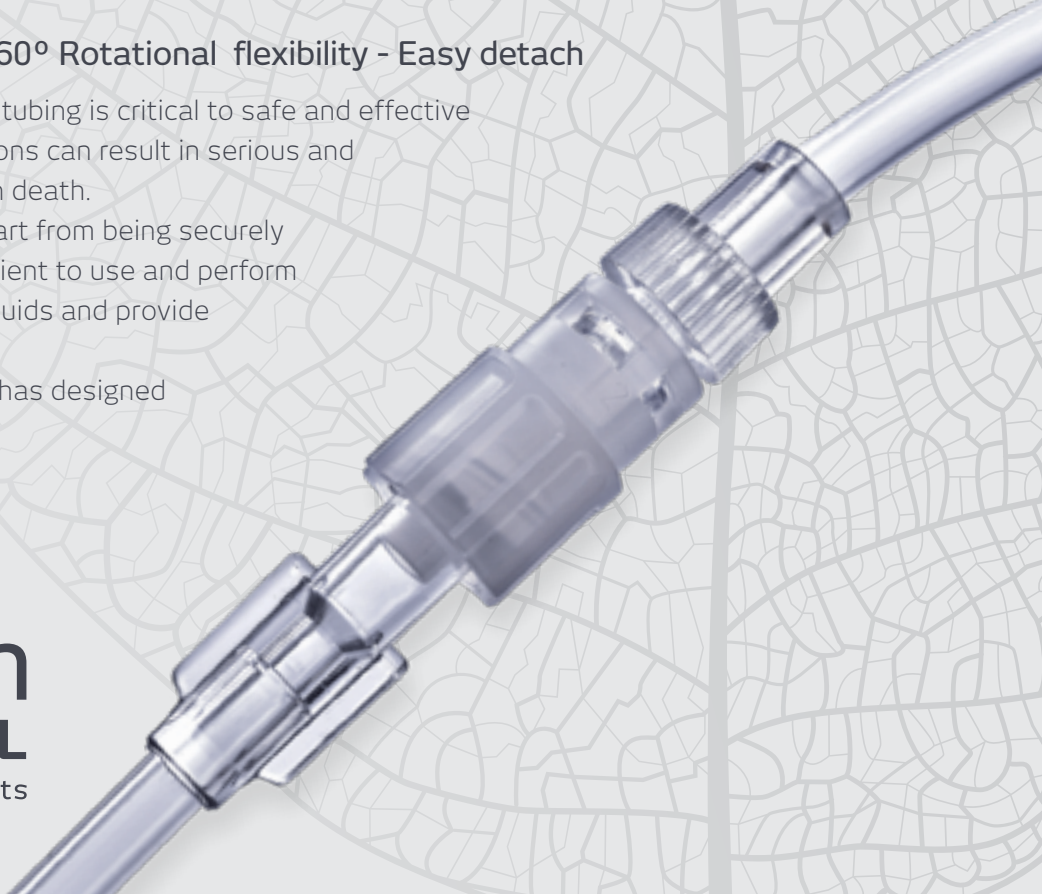

IV Therapy & Anesthesia



 Elcam Devices  Elcam Stopcocks  Elcam Components

Safe2 Rotator™ Luer Connector

The connection you can rely on

The Safe2 Rotator™ spinning lock feature provides ease of use and an extra layer of passive protection to your most critical fluid managing applications, without altering the traditional usability and common practice:

- Prevents accidental disconnections
- Facilitates convenient approach to injection and sampling ports (for example stopcock's side port and Y-site) during fluid administration
- Enables rotation of coupled luer lock connections and stopcocks
- Reduces tube kinks and twists in the tubing set (Kinks in the tubing can reduce the flow rates or cause blood hemolysis during transfusion or sampling)
- Prevents leakage due to over-tightening



Especially convenient in use with needleless connectors providing rotational flexibility of the distal part of the set.

Elcam Medical's Safe2 Rotator™ – secured connection along the entire treatment

The Safe2 Rotator™ Luer Connector is available in both stand-alone and stopcock configurations. Stopcock configuration is CE approved for non-sterile products and 510(k) cleared



www.elcam-medical.com

Head Office/Manufacturer
Elcam Medical

BarAm 1386000,
Israel
Tel: 972 4 698 8120/1/2
Fax: 972 4 698 0777
sales@elcam.co.il

USA
Elcam Medical Inc.

2 University Plaza, Suite 620
Hackensack, NJ 07601
Tel: +1 201 4571120
Fax: +1 201 4571125
info@elcam-medical.com

Europe
Elcam Medical Europe

Via Emilia Romagna 15,
41012 Carpi (Modena) Italy
Tel: +39 059 697710
Fax: +39 059 697705
sales.eu@elcam-medical.com

China
Elcam China Office

Room 3309, Jinhai Mansion
No. 89 of Anshun Road,
Shanghai, PRC. Zip: 200051
Tel: +86 21 62829967
Fax: +86 21 62829967
johnliu02@126.com